

Jesse M. Carter
P.O. Box 13275
Tampa, FL 33681-3275

1743

PLACE STICKER AT TOP OF ENVELOPE TO THE RIGHT
OF THE RETURN ADDRESS, FOLD AT DOTTED LINE
CERTIFIED MAIL™



7004 2890 0000 9149 3020


RECEIVED
MAR 07 2005
USPTO MAIL CENTER



U.S. POSTAGE
PAID
TAMPA, FL
33611-05
NO POSTAGE
NECESSARY
IF MAILED
IN THE
UNITED STATES
\$4.88
00067333-06

U.S. Patent and Trademark Office
Commissioner for Patents
Attn: AIA Unit - Helena G. Galka
P.O. Box 1450
Alexandria, VA 22313-1450

1743

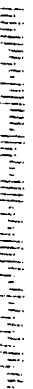

3-7-5

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

In RE:

APPLICANT(S): JESSE M. CARTER
ROBERT L. IMMEKUS

GROUP ART UNIT: 1743

S.N.: 10/051,845

EXAMINER: GAKH, YELENA G.

FILED: 01/17/2002

FOR: METHOD OF DETECTING OXIDIZING ADULTERANTS IN URINE

Commissioner of Patents and Trademarks
BOX NON FEE AMENDMENT
Washington, D.C. 20231

Madam:

In response to the Office Action dated 2-15-05 applicant respectfully requests reconsideration of rejection of claims 1-10 based on the following amendments and remarks.


Amendments and Remarks

Responses to detailed action page 2 of the Office Action dated 2-15-05.

1. Applicant begs your pardon as regards the non-compliance of the amendments sent to you 6-1-04. Applicant respectfully requests replacement of paragraph #1 on page 6, paragraph #2 on page 8, paragraph #1 on page 12, and paragraph 2 on page 13 according to 37 CFR 1.121. A copy of each paragraph showing corrections and a copy of each replacement paragraph are enclosed.

4. Applicant begs your pardon as regards the non-compliance of the claims sent to you 6-1-04. Applicant respectfully submits a complete list of claims 1-10, with corrections illustrated and appropriate status identifiers for each claim.

Sincerely,

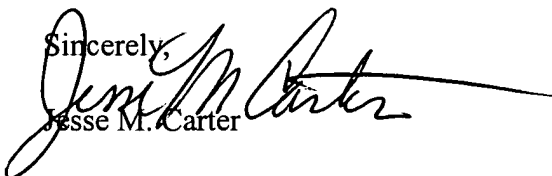
Jesse M. Carter

Corrected copy of page 6 paragraph #1 showing deletions and additions.

More precisely described, the invention is an automated method for detecting the presence of oxidants in a urine sample comprising placing an aliquot of the urine in a first automated analyzer sample cup, placing a standard of known concentration of oxidant in a second automated analyzer sample cup, placing the cups in a sample tray within the automated analyzer, transferring the urine from the first sample cup and the standard from the second sample cup to discrete cuvettes mounted within the automated analyzer, injecting an aqueous reagent composition comprising an acid at 0.01 N concentration or higher and one or more phenylamine chromogenic indicators from the group consisting of

N,N, N,N, N',N'-tetramethyl-1,4-phenylenediamine, N,N-diethyl-1,4-phenylenediamine, 2,3,5,6-tetramethyl-1,4-phenylenediamine, N,N-dimethyl-p-phenylenediamine, 2,4,6-trimethyl-1,3-phenylenediamine, N,N, N,N, N',N'-tetramethylbenzidine, 3,3,5,5-tetramethylbenzidine, N,N, N,N, N',N'-tetramethyl-4,4-diaminestilbene and O-tolidine.

into cuvettes and mixing and determining the sample-reagent mixture's absorbance with the automated analyzer's spectrophotometer at a preprogrammed wavelength between 400 and 700 nm at a preprogrammed time between 12 seconds and 600 seconds after the phenylamine chromogenic indicator is added to the cuvettes and comparing the absorbance of the urine-reagent mixture to the absorbance of the standard-reagent mixture and thereby determining if the sample has an abnormal amount of oxidant compound present. The reagent may also contain potassium iodide to enhance sensitivity to bleach and iodine containing oxidants. This enhanced sensitivity is important because bleach and some other halides will bind to nitrogen containing compounds in urine.

Corrected copy of page 8 paragraph #2 showing deletions and additions.

More precisely described, the invention is an automated method for detecting the presence of oxidants in a urine sample comprising placing an aliquot of the urine in a first automated analyzer sample cup, placing a standard of known concentration of oxidant in a second automated analyzer sample cup, placing the cups in a sample tray within the automated analyzer, transferring the urine from the first sample cup and the standard from the second sample cup to discrete cuvettes mounted within the automated analyzer, injecting a first aqueous reagent composition comprising potassium iodide and one or more buffering compounds into the cuvettes, injecting a second aqueous reagent composition comprising an acid at 0.01 N concentration or higher and one or more phenylamine chromogenic indicators from the group consisting of

N,N, N,N, N',N'-tetramethyl-1,4-phenylenediamine,

N,N-diethyl-1,4,-phenylenediamine, 2,3,5,6-tetramethyl-1,4-phenylenediamine,

N,N-dimethyl-p-phenylenediamine, 2,4,6-trimethyl-1,3-phenylenediamine,

N,N, N,N, N',N'-tetramethylbenzidine, 3,3,5,5-tetramethylbenzidine,

N,N,N,N,N',N'-tetramethyl-4,4-diaminestilbene and O-tolidine into the cuvettes and mixing and determining the sample-reagent mixture's absorbance with the automated analyzer's spectrophotometer at a preprogrammed wavelength between 400 and 700 nm at a preprogrammed time between 12 seconds and 600 seconds after the phenylamine chromogenic indicator is added to the cuvettes and, comparing the absorbance of the urine-reagent mixture to the absorbance of the standard-reagent mixture and thereby determining if the sample has an abnormal amount of oxidant compound present.

Corrected copy of page 12 paragraph #1 showing deletions and additions.

EXAMPLE IV

Prepare a solution containing:

R1

- a) 11.75 g Potassium Iodide
- b) 34.0 g Sodium Acetate
- c) 2.94 mLs 5.0 N Sodium Hydroxide
- d) QS to 1 liter with water.

R2

- a) 0.1 g DEPD (Diethyl-1,4-phenylenediamine sulfate)
- b) 0.333 g N,N, ~~N,N~~, N',N'-Tetramethyl-1,4-phenylenediaminedihydrochloride
- c) 6.9 mLs Phosphoric Acid
- d) QS to 1 liter with water.

This formulation is added to samples at a ratio of 1 to 7 to 7 (e.g. 18 μ L to 130 μ L to 130 μ L). This assay would be calibrated with 150 mg/L of nitrite standard and absorbance measured at 570 nm. This formula has good sensitivity to bleach, nitrite, chromate, iodate/iodic acid, and peroxide/peroxidase. The potassium iodide acts to intensify the reactivity of bleach and iodic acid.

Corrected copy of page 13 paragraph #2 showing deletions and additions.

EXAMPLE V

Prepare a liter of solution containing:

A. 0.25 g N,N, N,N, N',N'-tetramethylbenzidine

B. 50 mLs 5 N Hydrochloric Acid

C. QS with water to make 1 liter.

This formulation is added to samples at a ratio of 13 to 1 (e.g., 15 μ L sample to 200 μ L reagent).

This assay would be calibrated with nitrite as the standard (200 mg/mL Nitrite) and absorbance measured at 415 nm. This formulation has good sensitivity to nitrites, chromate, and peroxide/peroxidase, but not to low levels of bleach and iodic acid. One could include 1.0 mL of Brij-30 per liter of reagent.

Replacement Paragraph #1, Page 6.

More precisely described, the invention is an automated method for detecting the presence of oxidants in a urine sample comprising placing an aliquot of the urine in a first automated analyzer sample cup, placing a standard of known concentration of oxidant in a second automated analyzer sample cup, placing the cups in a sample tray within the automated analyzer, transferring the urine from the first sample cup and the standard from the second sample cup to discrete cuvettes mounted within the automated analyzer, injecting an aqueous reagent composition comprising an acid at 0.01 N concentration or higher and one or more phenylamine chromogenic indicators from the group consisting of

N,N,N',N'-tetramethyl-1,4-phenylenediamine, N,N-diethyl-1,4-phenylenediamine,

2,3,5,6-tetramethyl-1,4-phenylenediamine, N,N-dimethyl-p-phenylenediamine,

2,4,6-trimethyl-1,3-phenylenediamine, N,N,N',N'-tetramethylbenzidine,

3,3,5,5-tetramethylbenzidine, N,N,N',N'-tetramethyl-4,4-diaminestilbene and O-tolidine.

into cuvettes and mixing and determining the sample-reagent mixture's absorbance with the automated analyzer's spectrophotometer at a preprogrammed wavelength between 400 and 700 nm at a preprogrammed time between 12 seconds and 600 seconds after the phenylamine chromogenic indicator is added to the cuvettes and comparing the absorbance of the urine-reagent mixture to the absorbance of the standard-reagent mixture and thereby determining if the sample has an abnormal amount of oxidant compound present. The reagent may also contain potassium iodide to enhance sensitivity to bleach and iodine containing oxidants. This enhanced sensitivity is important because bleach and some other halides will bind to nitrogen containing compounds in urine.

Replacement Paragraph #2, Page 8.

More precisely described, the invention is an automated method for detecting the presence of oxidants in a urine sample comprising placing an aliquot of the urine in a first automated analyzer sample cup, placing a standard of known concentration of oxidant in a second automated analyzer sample cup, placing the cups in a sample tray within the automated analyzer, transferring the urine from the first sample cup and the standard from the second sample cup to discrete cuvettes mounted within the automated analyzer, injecting a first aqueous reagent composition comprising potassium iodide and one or more buffering compounds into the cuvettes, injecting a second aqueous reagent composition comprising an acid at 0.01 N concentration or higher and one or more phenylamine chromogenic indicators from the group consisting of N,N,N',N'-tetramethyl-1,4-phenylenediamine,

N,N-diethyl-1,4-phenylenediamine, 2,3,5,6-tetramethyl-1,4-phenylenediamine,

N,N-dimethyl-p-phenylenediamine, 2,4,6-trimethyl-1,3-phenylenediamine,

N,N,N',N'-tetramethylbenzidine, 3,3',5,5'-tetramethylbenzidine,

N,N,N',N'-tetramethyl-4,4'-diaminostilbene and O-tolidine into the cuvettes and mixing and determining the sample-reagent mixture's absorbance with the automated analyzer's spectrophotometer at a preprogrammed wavelength between 400 and 700 nm at a preprogrammed time between 12 seconds and 600 seconds after the phenylamine chromogenic indicator is added to the cuvettes and, comparing the absorbance of the urine-reagent mixture to the absorbance of the standard-reagent mixture and thereby determining if the sample has an abnormal amount of oxidant compound present.

EXAMPLE IV

Prepare a solution containing:

R1

- a) 11.75 g Potassium Iodide
- b) 34.0 g Sodium Acetate
- c) 2.94 mLs 5.0 N Sodium Hydroxide
- d) QS to 1 liter with water.

R2

- a) 0.1 g DEPD (Diethyl-1,4-phenylenediamine sulfate)
- b) 0.333 g N,N,N',N'-Tetramethyl-1,4-phenylenediaminedihydrochloride
- c) 6.9 mLs Phosphoric Acid
- d) QS to 1 liter with water.

This formulation is added to samples at a ratio of 1 to 7 to 7 (e.g. 18 μ L to 130 μ L to 130 μ L). This assay would be calibrated with 150 mg/L of nitrite standard and absorbance measured at 570 nm. This formula has good sensitivity to bleach, nitrite, chromate, iodate/iodic acid, and peroxide/peroxidase. The potassium iodide acts to intensify the reactivity of bleach and iodic acid.

EXAMPLE V

Prepare a liter of solution containing:

A. 0.25 g N,N,N',N'-tetramethylbenzidine

B. 50 mLs 5 N Hydrochloric Acid

C. QS with water to make 1 liter.

This formulation is added to samples at a ratio of 13 to 1 (e.g., 15 μ L sample to 200 μ L reagent).

This assay would be calibrated with nitrite as the standard (200 mg/mL Nitrite) and absorbance measured at 415 nm. This formulation has good sensitivity to nitrites, chromate, and peroxide/peroxidase, but not to low levels of bleach and iodic acid. One could include 1.0 mL of Brij-30 per liter of reagent.